





Single Infusion Line



2 - Way TIVA Set

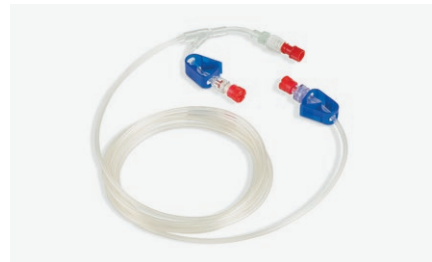

2 - Way Infusion Set

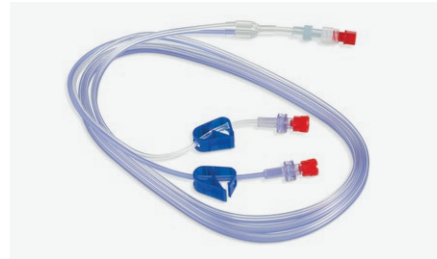

Dual TIVA Set, non-DEHP



3 - Way TIVA Set



4 - Way Triple Plus TIVA Set


Ordering information

Product Code	Product Description	Length	Box Quantity
6611	Single Infusion Line	2.0m	50
6616	Single Infusion Line	6.0m	25
6202	2 - Way Infusion Set	2.5m	50
6610	Dual TIVA set	2.0m	50
6604	2 - Way TIVA Set	2.0m	50
6603	3 - Way TIVA Set	2.0m	50
6605	3 - Way TIVA Set	2.5m	25
6606	3 - Way TIVA Set	3.0m	25
6607	3 - Way TIVA Set	6.0m	10
6630	4 - Way Triple Plus TIVA Set	2.2m	25
6631	4 - Way Triple Plus TIVA Set	3.0m	25
6633	Triple MRI Extension Set	6.0m	10


Clinical data

Guidelines for the safe practice of total intravenous anaesthesia (TIVA) Joint Guidelines from the Association of Anaesthetists and the Society for Intravenous Anaesthesia A. F. Nimmo, A. R. Absalom, O. Bagshaw, A. Biswas, T. M. Cook A. Costello, S. Grimes, D. Mulvey, S. Shinde, T. Whitehouse and M. D. Wiles.

Tips and Tricks to optimize Total Intravenous Anesthesia L. Barvais, F. A. Lobo, F. H. M. Engbers, M. G. Irwin, T. W. Schnider and S. Schraag.

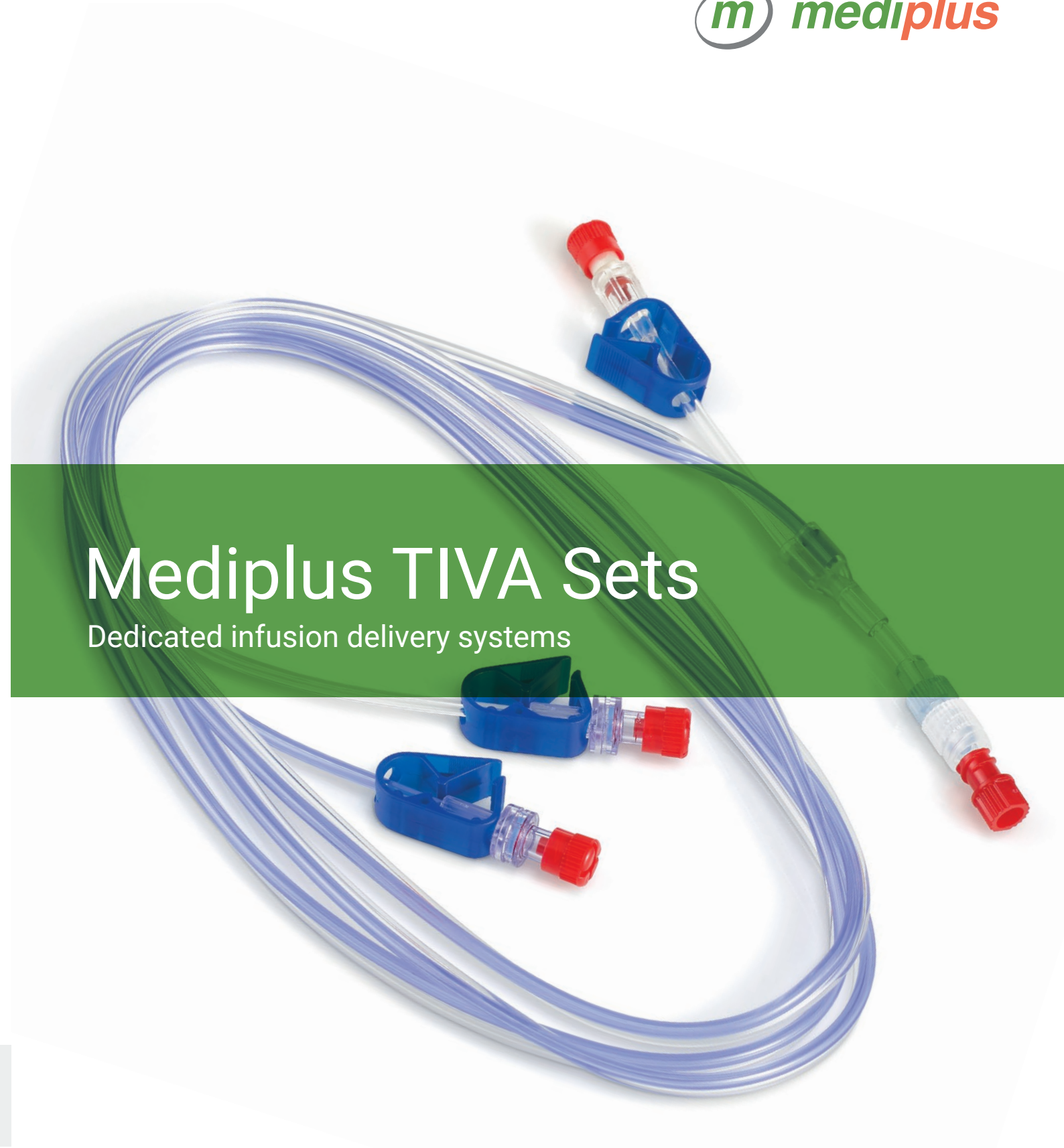

Mediplus are experts in the delivery of high performance medical devices that support positive patient outcomes in the following fields -

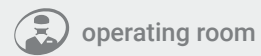

Mediplus Ltd.
 +44 (0) 1494 551200
 help@mediplusuk.com
 The Gateway Centre, Coronation Road,
 High Wycombe HP12 3SU
 www.mediplusuk.com



INVESTORS IN PEOPLE™
 We invest in people Silver



Mediplus TIVA Sets

Dedicated infusion delivery systems

Mediplus TIVA Sets were designed in conjunction with anaesthetists to meet the highest TIVA standards as outlined in the SIVA, Euro SIVA and AAGBI guidelines. These high-performance sets offer a wide range of benefits including –

Enhanced performance

- Using a dedicated line reduces the risk of leakage associated with multiple connections
- Integrated anti-siphon valve and anti-reflux valve prevent free flow and back flow
- Wide bore IV gravity line enables high fluid flow to ensure continuous infusion
- Minimal dead space designed to help prevent potential overdose in cases of accidental bolus
- Micro bore tubing to prevent kinking and occlusion
- Siamese tubing on double and triple infusion lines to prevent entanglement

Greater patient comfort

- The soft distal end minimises stress at the cannula site to help avoid displacement

Cost savings

- Micro bore tubing minimises post-operative drug wastage

Wide range

- Single, double and triple infusion lines available in a range of lengths

“We used the Mediplus TIVA giving set for over 15 years without incident or failure. As a cost saving measure we changed to an alternative giving set, but within four months we found that product to be inferior and we have now returned to our tried and trusted product”

Bruce Rodda, Consultant Anaesthetist
Mid Essex Hospital Services NHS Trust


Key features



Wide bore gravity line to allow high flow of IV fluid



Siamese configuration to prevent entanglement of tubing

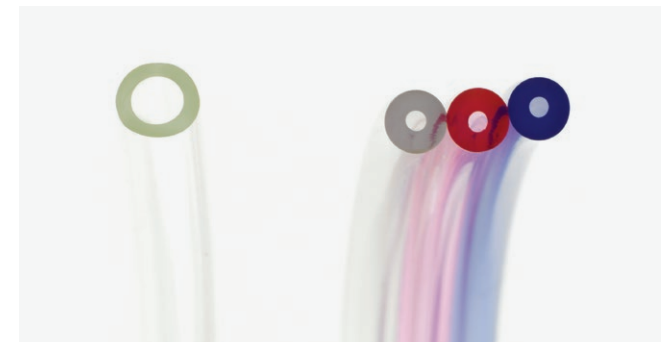

Micro bore siamese tubing for a low priming volume but large external diameter to prevent kinking/occlusion



Anti-Siphon valves to prevent back flow and free-flow



Anti-reflux valve to prevent back flow



Soft distal end to provide minimal stress at the cannula site